



# PT400 PHYSICAL THERAPY TABLE

Rock solid support for better results




Shown in Pewter fabric



Chair  
Position



Powered Center Break  
0° to 45°

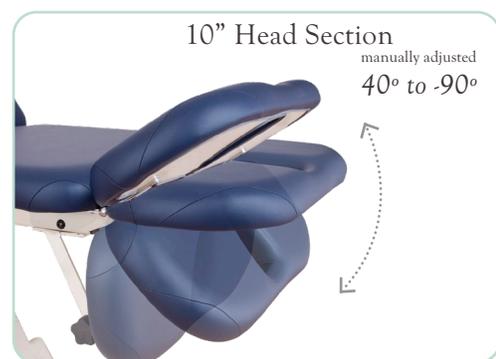

10" Head Section  
manually adjusted  
40° to -90°

## Features

- 550 lb weight capacity
- Powered center break
- Can be used as a chair or table
- Uniform stability maintained in 16"-34"  
electric height range



Twin beam design lifts vertically, taking up much less space than cantilevered designs.

Visit [www.oakworksmed.com](http://www.oakworksmed.com) for a quote!

# PT400

#70508-T\* 29" W 1.75" Firm Response    #73393-T\* 29" W 2.5" Comfort  
 #71775-T\* 31" W 1.75" Firm Response    #73392-T\* 31" W 2.5" Comfort

FEATURES	
3 Section Top	29" (74cm) or 31" (79cm) wide x 73" (185cm) long
Backrest	39" (99cm) long
Seat	23" (58cm) long
Head	11" (28cm) long
Motion 1	Electric Height Range: 16" - 34" (40-86cm) (1.75" foam)
Motion 2	Electric Center Break 0° - 45°
Motion 3	Manual Lift-Assist Backrest 0° - 85°
Motion 4	Manual position head section 40° to -90°
Controls	Multi-function foot control
Weight Capacity	550 lb (250kg)
Product Weight	198 lb (90 kg)
Padding	1.75" (4cm) Firm Response or 2.5" (6cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Locking Twin Casters (one end of the table)	To easily move your table in the room
Side Rails	Drop away (available on 31" only)
Second Foot Control	Controls height & center break
Hand Control	Controls height & center break

ACCESSORIES <i>order any time.</i>	
Universal Paper Roll Holder	Mounts to installed accessory mounting plate
Face Hole Filler	Fits securely in cut out



Optional Locking  
Twin Casters

## STABLE TWIN BEAM CONSTRUCTION

Stiffer than cantilever designs and with fewer moving parts, our unique design lifts vertically, taking up less room than a cantilevered design.



## FIRM RESPONSE FOAM

Firm Response foam delivers a firm yet comfortable surface that maximizes ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design.



## LIFT ASSIST BACKREST

Lift assist backrest features two powerful pneumatic cylinders that assist with adjusting the backrest up or down.



## TERRATOUGH™ FABRIC

Coal-T01	White-T26	Heron-T05	Earth-T07
Opal-T12	Stone-T13	Sapphire-T17	Clay-T18
Ocean-T20	Pewter-T29	Sage-T22	Blue Grass-T24
Orchid-T19	Espresso-T27	Sky Blue-T04	Ruby-T08

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

717.759.3146 [www.oakworksmid.com](http://www.oakworksmid.com) email: [medsales@oakworksmid.com](mailto:medsales@oakworksmid.com)